


Overbath Stretcher Pre-Installation Requirements

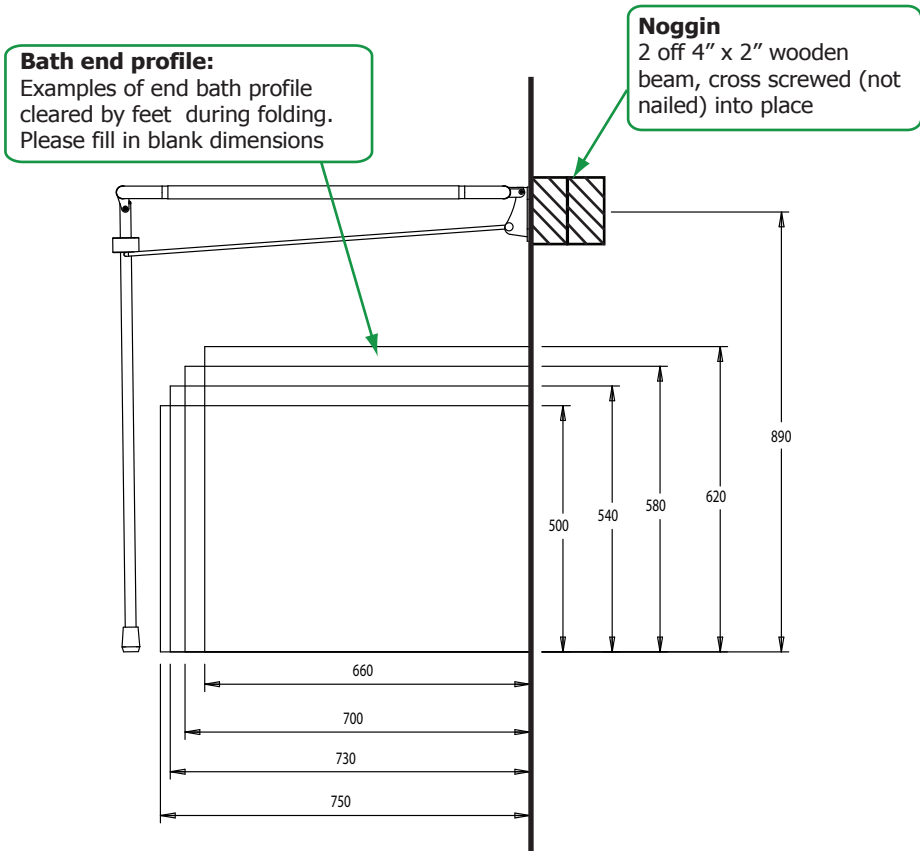

Wall Specification

The changing table should be fixed to a load bearing wall (ie. strong enough to carry a 100kg load on wall fixing) which is flat and free from obstruction from 700mm to 1800mm above the floor. If the changing table must be mounted on a stud wall, the wall will need special reinforcement in the form of either noggins placed in the wall as shown above, or 25mm (1 inch) marine or exterior plywood used to cover the entire fixing area.

We understand that any waiting time will be charged to us and should the site require a second visit due to the above requirements not being met an additional £350.00 per visit will be charged and we agree to pay those charges.

Signed Date

Print name Job number

Site address